oserio®

www.oserio.com

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註：體脂肪%的標準參考文獻：Daniel P. Williams, Scott B. Going, Timothy G. Lohman, David W. Harsha, Satharnur R Srinivasan, Larry S. Webber, Gerald S. Berenson, Body fatness and risk for elevated blood pressure, total cholesterol, and serum lipoprotein ratios in children and adolescents American journal of public health, March 1992, Vol. 82, No.3

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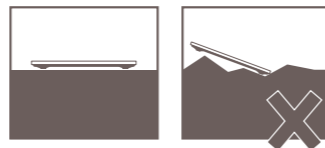
Guidance and manufacturer's declaration-electromagnetic emissions

The FFP-330A; FFP-330B; FFP-330C; FFP-330D; FFP-331 is intended for use in the electromagnetic environment specified below. The customer or the user of the FFP-330A; FFP-330B; FFP-330C; FFP-330D; FFP-331 should assure that it is used in such an environment.

Emission test	Compliance	Electromagnetic environment-guidance
RF emissions CISPR 11	Group 1	The FLG-351; FEG-113; FEG-123 uses RF energy only for its internal function. Therefore, its RF emissions are very low and are not likely to cause any interference in nearby electronic equipment.
RF emissions CISPR 11	Class B	The FLG-351; FEG-113; FEG-123 is suitable for use in all establishments, including domestic establishments and those directly connected to the public low-voltage power supply network that supplies buildings used for domestic purposes.
Harmonic emissions IEC 61000-3-2	Not applicable	
Voltage fluctuations/flicker emissions IEC 61000-3-3	Not applicable	

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Guidance and manufacturer's declaration-electromagnetic immunity			
The FLG-351; FEG-113; FEG-123 is intended for use in the electromagnetic environment specified below. The customer or the user of the FLG-351; FEG-113; FEG-123 should assure that it is used in such an environment.			
Immunity test	IEC 60601 test level	Compliance level	Electromagnetic environment-guidance
Electrostatic discharge(ESD) IEC 61000-4-2	± 6 kV contact ± 8 kV air	± 6 kV contact ± 8 kV air	Floors should be wood, concrete or ceramic tile. If floors are covered with synthetic material, the relative humidity should be at least 30%.
Electrical fast transient/burst IEC 61000-4-4	± 2kV for power supply lines + 1kV for input/output lines	± 2kV for power supply lines Not applicable	Mains power quality should be that of a typical commercial or hospital environment.
Surge IEC 61000-4-5	± 1kV line(s) to line(s) ± 2kV line(s) to earth	± 1kV differential mode Not applicable	Mains power quality should be that of a typical commercial or hospital environment.

Guidance and manufacturer's declaration-electromagnetic immunity			
Immunity test	IEC 60601 test level	Compliance level	Electromagnetic environment-guidance
Voltage Dips, short interruptions and voltage variations on power supply input lines IEC 61000-4-11	<5% UT(>95% dip in UT) for 0.5 cycle 40% UT(60% dip in UT) for 5 cycles 70% UT(30% dip in UT) for 25 cycles <5% UT(>95% dip in UT) for 5 s	Not applicable Not applicable Not applicable Not applicable	Mains power quality should be that of a typical commercial or hospital environment. If the user of the FLG-351; FEG-113; FEG-123 requires continued operation during power mains interruptions, it is recommended that the FLG-351; FEG-113; FEG-123 be powered from an uninterruptible power supply or a battery.
Power frequency (50/60 Hz) magnetic field IEC 61000-4-8	3 A/m	3 A/m	The FLG-351; FEG-113; FEG-123 power frequency magnetic fields should be at levels characteristic of a typical location in a typical commercial or hospital environment.
NOTE: UT is the a.c. mains voltage prior to application of the test level.			

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Guidance and manufacturer's			declaration-electromagnetic immunity	
Immunity test	IEC 60601 test level	Compliance level	Electromagnetic environment-guidance	
Conducted RF IEC 61000-4-6	3 Vrms 150 KHz to 80 MHz	Not applicable	<p>Portable and mobile RF communications equipment should be used no closer to any part of the FLG-351; FEG-113; FEG-123 including cables, than the recommended separation distance calculated from the equation applicable to the frequency of the transmitter.</p>	
Radiated RF IEC 61000-4-3	3 V/m 80MHz to 2.5 GHz	3 V/m	<p>Recommended separation distance: $d = 1.2 \sqrt{P}$ $d = 1.2 \sqrt{P}$ 80MHz to 800 MHz $d = 2.3 \sqrt{P}$ 800MHz to 2.5 GHz Where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer and d is the recommended separation distance in metres (m). Field strengths from fixed RF transmitters, as determined by an electromagnetic site survey^a, should be less than the compliance level in each frequency range^b.</p> <p>Interference may occur in the vicinity of equipment marked with the following symbol: </p>	

NOTE1:

At 80 MHz and 800 MHz, the higher frequency range applies.

NOTE2:

These guidelines may not apply in all situations.

Electromagnetic propagation is affected by absorption and reflection from **structures, objects and people**.

- A. Field strengths from fixed transmitters, such as base stations for radio (cellular/cordless) telephones and land mobile radios, amateur radio, AM and FM radio broadcast and TV broadcast cannot be predicted theoretically with accuracy. To assess the electromagnetic environment due to fixed RF transmitters, an electromagnetic site survey should be considered. If the measured field strength in the location in which the **FLG-351; FEG-113; FEG-123** should be observed to verify normal operation. If abnormal performance is observed, additional measures may be necessary, such as re-orienting or relocating the **FLG-351; FEG-113; FEG-123**. Over the frequency range 150 kHz to 80 MHz, field strengths should be less than 3
- B. V/m.

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Recommended separation distance between portable and mobile RF communications				equipment and the FLG-351; FEG-113; FEG-123	
The FLG-351; FEG-113; FEG-123 is intended for use in an electromagnetic environment in which radiated RF disturbances are controlled. The customer or the user of the FLG-351; FEG-113; FEG-123 can help prevent electromagnetic interference by maintaining a minimum distance between portable and mobile RF communications equipment (transmitters) and the FLG-351; FEG-113; FEG-123 as recommended below, according to the maximum output power of the communications equipment.				For transmitters rated at a maximum output power not listed above, the recommended separation distance d in meters (m) can be estimated using the equation applicable to the frequency of the transmitter, where p is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer.	
Rated maximum output power of transmitter (W)	Separation distance according to frequency of transmitter (m)			NOTE1: At 80 MHz and 800 MHz, the separation distance for the higher frequency range applies. NOTE2: These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.	
	150 kHz to 80 MHz $d = 1.2 \sqrt{P}$	80 MHz to 800 MHz $d = 1.2 \sqrt{P}$	800 MHz to 2.5 GHz $d = 2.3 \sqrt{P}$		
0.01	N/A	0.12	0.23		
0.1	N/A	0.38	0.73		
1	N/A	1.2	2.3		
10	N/A	3.8	7.3		
100	N/A	12	23		

1

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
Explanation of the graphic symbols:

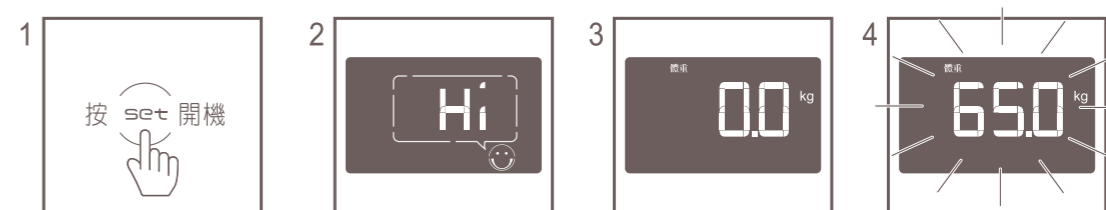SN-135-00001	Designation of the serial number of every device, applied at the device. (Number as an example)		Carefully read this operation manual before setup and commissioning, even if you are already familiar with oserio scales.
	"Please note the accompanying documents" or "Observe operating instructions"		Dispose of old appliances separately from your household waste! Instead, take them to communal collection points.
 Chander Electronic Co., Ltd. No.103, Guozhong Rd., Dali Dist., Taichung City, 412 Taiwan (R.O.C.)	Identification of manufacturer of medical product including address		Transport and storage temperature limit indicating the upper and the lower limit. (Transport and storage temperature on packaging)
	"Electro-medical appliance" with attachment for type BF		

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